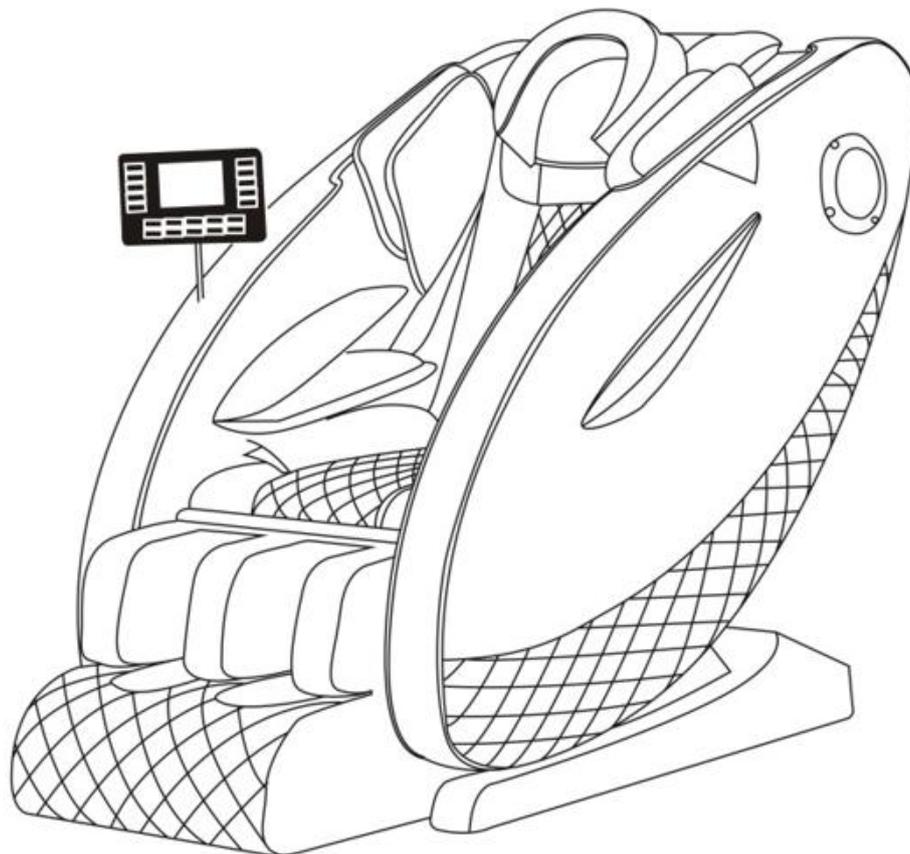
BODY FLEX[®]

MAS 6 MASSAGE CHAIR

art. BF-MASCH6

INSTRUCTION MANUAL

- The whole machine has beautiful lines,bright colors, humanized design and green environmental protection materials.
- This product is designed with 2 groups of movements and 8 massage heads, which can knead and massage the neck,back,waist of human body.
- The whole machine is equipped with about 20 air bags divided into legs,arms, cushions and shoulders to freely display the function of air bag finger pressing and massage.
- It can manually adjust the massage movement of neck,back,waist.Various massage techniques can be changed freely,and three-level adjustment can be realized.
- The recline angle of the back can be adjusted according to the needs,and the free adjustment of 115°to 155°can be realized.When the angle of the back is adjusted to the minimum,the zero gravity function can be realized.
- One key heating,constant temperature 35°to 45°C.
- The massage time can be adjusted according to the needs,with 5 minutes plus and minus steps,and a massage journey can be set for 15 minutes,20 minutes, 25minutes and 30 minutes.
- ▶ Specially designed sole rollers.

ATTENTION: The description on the website of the seller is the leading and valid reference. The manufacturer reserves the right to make changes to the specifications. Some of the described features may not be available on all chair models.

2.Names of components

Usage

1.Installation Method of Connection Line

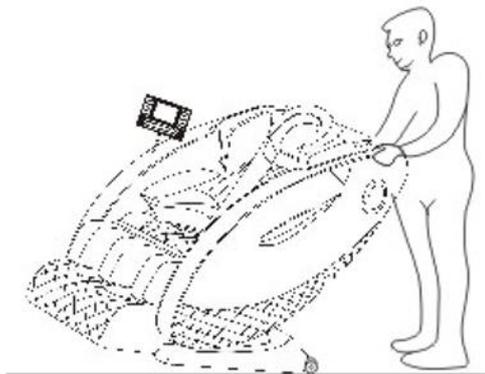
Figure 1

Connect the type plug of the power cord to the plug of the electrical box(as shown in Figure 1) and the other end to the power socket.

Movement method

When moving to support the lower part of the pad, lift the front move. **Attention!**
Do not move people
Do not move with small pulleys on floors that are prone to damage

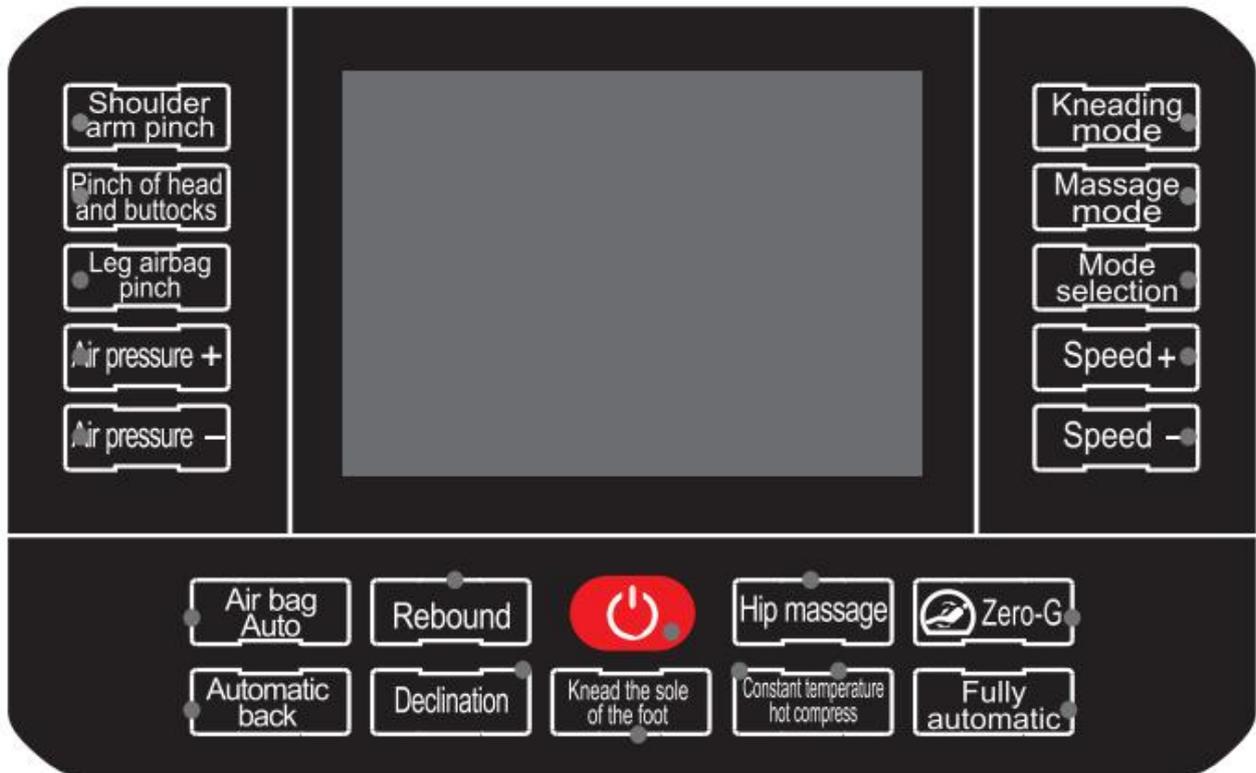
Placement of massage chair

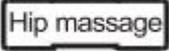
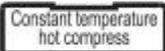
Please leave the proper space around the massage chair when preparing for the massage.

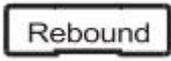
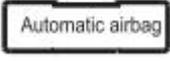
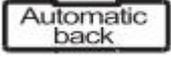
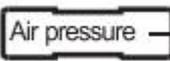
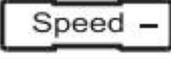
- (1)The place where the massage chair is placed should have enough space for the chair to recline and the foot stand to rise.
- (2)Make sure there is at least 80cm of space in the back of the massage chair and at least 50cm in the front
- (3)Make sure there are no obstacles in front and behind the massage chair
- (4)In order to prevent the ground from being damaged, please place a carpet under the massage chair.so as not to damage the ground

3.Manual operation instructions

Remote Control Keys

	<p>Press one time start/Timing/Long press to shut down Press one time start: When the power is on, press this button and the massage chair enters the automatic mode. The default working time is 15 minutes. Timing: When the massage chair is working, press this button, the massage chair enters the timed mode, and each time is increased by 5 minutes. Long press to shut down: When the massage chair is working, press and hold this button until the display screen is turned off.</p>
	<p>Full automatic: press this key to set five modes, and select F1/F2/F3/F4/F5.</p>
	<p>Zero gravity: press this key to recline / reset the massage chair.</p>
	<p>Hip massage: press this key to turn on / off the cushion massage function.</p>
	<p>Constant temperature hot compress: press this key to turn on /off the constant temperature hot compress function of calf and waist.</p>
	<p>Foot kneading: press this key to turn the foot kneading function on / off.</p>

	Pick up: long press this key to adjust the reset and pick up function of massage chair (long press is effective).
	Recline: long press this key to adjust the recline function of massage chair (long press is effective).
	Automatic airbag: press this key to turn on / off the whole body airbag squeeze massage.
	Back automatic: press this key to turn back movement massage on / off.
	Shoulder arm pinch: press this key to turn on /off shoulder arm airbag pinch massage.
	Head and arm pinching: press this key to turn on / off the head and hip airbag pinching massage.
	Leg airbag pinch: press this key to turn on / off leg airbag pinch massage.
	Air pressure +: the air bag extrusion force increases.
	Air pressure -: the pressure of the air bag is reduced.
	Kneading mode: press this key to turn on / off neck kneading massage.
	Massage mode: press this key to turn on / off waist massage.
	Mode selection: press this key to select the direction of the massage movement (positive/negative).
	Speed +: the speed of the massage movement increases.
	Speed -: the speed of the massage movement slows down.

4.Maintenance and Maintenance

- 1.Please use neutral detergent and soft cloth for daily cleaning,and do not use corrosive detergent.
- 2.If the power cord is damaged,it must be replaced by the full-time personnel of the maintenance shop designated by the manufacturer to avoid danger.
- 3.When not using this product or cleaning,please unplug the power supply to avoid danger
- 4.Do not use any sharp objects to cut the product surface.
- 5.Do not put the product in a humid environment,but in a dry and ventilated place.
- 6.Do not press or press the product hard.

5.Matters needing attention

In order to ensure the safety and correct use of the massage chair,please refer to the operation manual in detail before using this product.

To reduce the risk of fire,burns,electric shock,or injury to others,follow these instructions:

- 1)Power supply other than AC 220V shall not be used.
- 2)When no one is using the machine,please be sure to remove the power plug of the massage chair from the socket.When not using the machine or before removing the device/parts,please also remove the power plug from the socket in time.
- 3)In order to avoid fire,electric shock or injury caused by overheating,please do not put the product under pillows,blankets and other covers.
- 4)If the product is used for children or the disabled,adult guardians are required to supervise and guide it.
- 5)Please be sure to use the product according to the usage described in the user manual. Do not use parts that are not recommended by the manufacturer.
- 6)Do not use this product in case of any of the following conditions:the power cord or plug is damaged,the operation is not normal,if it falls to the ground or other damages, or falls into water,please send the product back to the service center for inspection and maintenance.
- 7)Do not drag the power cord to move the product or use the power cord as a handle.
- 8)Keep a safe distance between the power cord and the heating surface.
- 9)Do not operate the product when the air outlet is blocked.Please keep the air outlet unobstructed to avoid being blocked by cotton chips,hair and other things.
- 10)Do not insert or drop objects into any openings of the product.
- 11)Do not use this product outdoors.
- 12)please do not operate the product in the place of spray product or oxygen filled place.
- 13)Cut off the power supply,please turn off all operation controls,and then unplug the power plug from the socket.
- 14)Do not cover the massage chair so as not to cause malfunction or fire due to excessive temperature.

6.Safety Protection

- Do not use the massage chair in excessively humid or dusty environment to avoid electric shock or mechanical failure.
- Do not use the massage chair when the room temperature exceeds 40°C.
- Do not place the massage chair near the electric heater, stove or direct sunlight.
- The massage chair is designed for personal use and is not suitable for any commercial use.

△ Careful

- Massage chair is only for indoor use.
- Use massage chair at low room temperature, do not suddenly adjust the temperature. It is recommended to gradually adjust the temperature to the normal range before use.
- If the massage chair is moved from a low temperature storage place to a warm environment for use, it may affect the normal operation of the machine due to the accumulation of moisture in the parts, or even cause parts failure. It is recommended that the product be placed at normal room temperature for 1 hour before use.

7. Warm Tips

- The power switch must be turned on before other functions can be selected.
- Children should use this massage chair under the supervision of adults to avoid equipment damage and children's injury due to improper operation.
- Please do not let children climb this massage chair to avoid children falling.
- Use the massage chair to adjust the strength of the chair massage, because the massage force is too small to play its due role in stimulating, too large prone to fatigue, and damage the skin, massage force and times from less to more, from light to heavy.
- Use the massage chair to relax the body and mind. In addition to thinking at home and abroad, you should be calm and relaxed.
- Massage chairs are not suitable for patients with heart Disease, hypertension, osteoporosis, local skin damage, ulcer, bleeding, tuberculosis, tumor, people who are too hungry, full, drunk or tired.
- In order to avoid the reduction of massage effect or injury, please do not use the massage chair together with other effective equipment such as electric blanket.
- Do not use the massage chair within 1 hour after eating to avoid adverse effects.
- Do not use the massage chair when your hands or other parts of your body are wet.
- Do not plug or unplug the power plug.
- Do not insert any articles or inflammables into the massage chair. This may cause electric shock, malfunction or fire.
- If you feel any discomfort during massage, please stop immediately.
- If the product is not used for self diagnosis or treatment and needs diagnosis or treatment, it shall be handled by professional medical personnel.
- Do not use this product if you are in medical condition and are pregnant or in physiological period.

Use the massage chair to master the massage time, 15 minutes for each massage, and no more than half an hour for each person (time) (when it reaches half an hour, let the massage chair rest for ten minutes before use). If it is used continuously for a long time, the overheat protection

device may be started and the power will be cut off automatically.

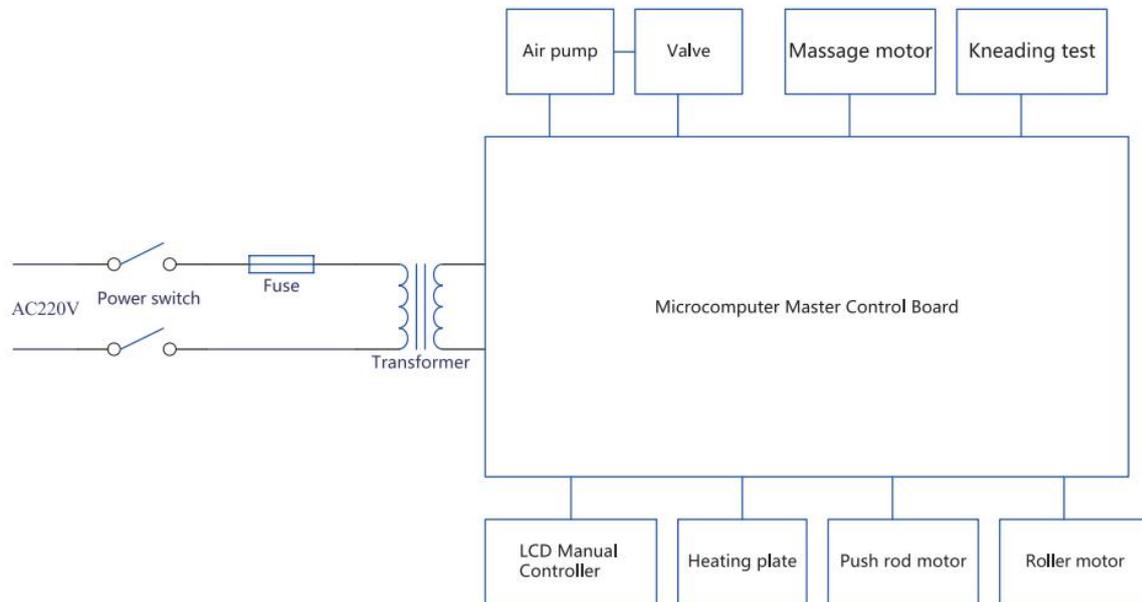
- The surface of the appliance is hot. People who are not sensitive to the heat must pay attention to it when using it.
- Do not directly massage elbow joint, knee, chest, abdomen, etc.

Not suitable for massage chair:

In case of any of the following conditions, please turn off the switch immediately, cut off the power supply and contact our company.

- accidentally spilling liquid on the massage chair may cause electric shock.
- the cloth cover of massage chair is damaged and the internal parts are exposed.
- in case of abnormal or severe pain during massage, please stop massage immediately and seek medical consultation.
- please stop using immediately if any mechanical failure or abnormal condition is detected during use.
- in case of current interruption in use, attention shall be paid to avoid user injury due to sudden power recovery.
- in case of lightning.

8.Circuit schematic diagram

9. Analysis and Removal of Common Faults

Failure phenomenon	Possible cause	Processing Method
Sound(working sound) *inflatable sound of airbags *bleed sound of airbags *changing-over sound on solenoid valve *air pump working sound	The sound produced by the product at work Sound of mechanical operation	This sound is normal No need to deal with it
The whole machine is not working at all	*The power plug is not inserted into the socket *Socket switch not open *The massage button is not pressed	Check the power plug and socket, ensure that the plug is in good contact with the socket, press the massage button.
The massage stopped suddenly in the process	This product has just completed 15 minutes of rated working time and automatically switched to standby mode.	If you would like to massage again, please do so again. Press the  key to turn on massage.
	Overheating protection for motors exceeding rated service time	Turn off, Cut the power, use the unit again after the motor is cool down

*The company reserves the right to improve products and instructions. If alteration is made without notice, the appearance of the product shall be in accordance with material object.

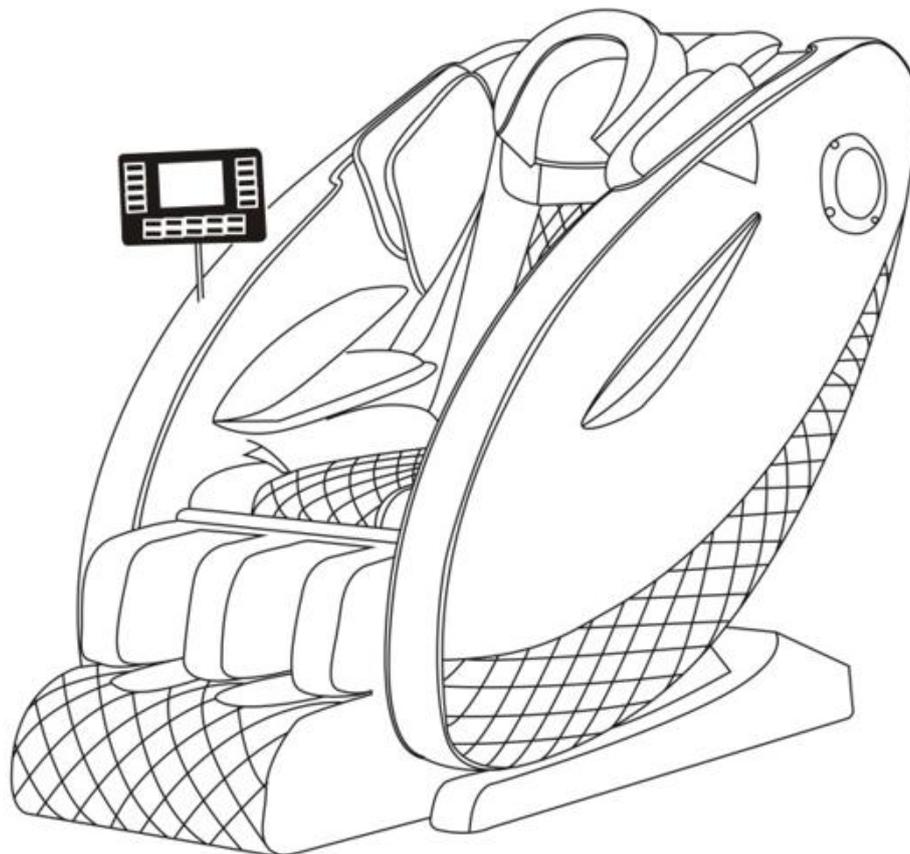
*This volume has been carefully checked, but it does not exclude the possibility of misprints or misunderstandings in the content. We are sorry to inform you that we reserve the right of interpretation.

Improper disposal of waste causes environmental pollution.

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